




PRESSALIT®

Product Item number	Nursing bench 2000, coated mesh cover, electrically height adjustable R8468000-93
Colour	<p>EAN:5708590317915 FDA:INQ/D124319</p> <p>White</p> 
Product description	<p>Discontinued per 31 Dec 2018 and replaced by R8478 Nursing bench 2000, powered, height adjustable 700 mm; from 300 mm -1000 mm, with wired hand control Foldable. Incl. water collection tray with integrated outlet. Outer measurements 880 mm x 1800 mm. Excl. mounting kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types.</p> <p>Discontinued per 31 Dec 2018 and replaced by R8478 Nursing bench 2000, powered, height adjustable 700 mm; from 300 mm -1000 mm, with wired hand control Foldable. Incl. water collection tray with integrated outlet. Outer measurements 880 mm x 1800 mm. Excl. mounting kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types.</p>
Adjustability Adjustment height Product capacity Speed upwards Noise level	<p>Height adjustable</p> <p>700 mm</p> <p>Maximum load: 175 kg. When the bench is loaded with 175 kg, each mounting bolt supports 23 kg <i>This product has been tested statically with a load of 440 kg as the IEC 60601-standard requires testing with safety factor 2.5.</i></p> <p>Speed up mm: 15 mm/s Speed down mm: 20 mm/s</p> <p>48 dB</p>

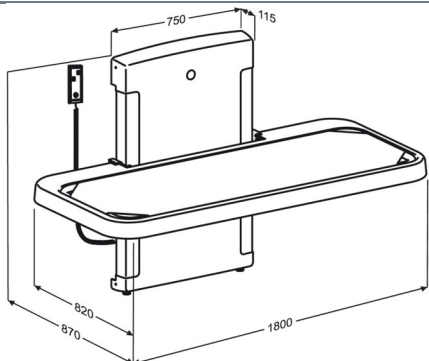

Product	Nursing bench 2000, coated mesh cover, electrically height adjustable
Item number	R8468000-93
Materials	The frame is made from powder-laquered stainless steel, the lying surface is a mesh cover. The water collection tray is made from polystyrene.
Ergonomy	<p>Pressalit Care offers a range of high quality, fixed or power adjusted nursing benches designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centres, institutions, hospitals and private homes.</p> <p>The nursing bench 3000 series provides for maximum safety and comfort. Besides an award winning design, the series includes comfortable and durable surfaces, including polyurethane foam and aluminum.</p> <p>The flexibility of the nursing bench provides for optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p>
Features	<p>A water collection tray with integrated water outlet prevents water from splashing on carers or the floor. The flexible hose attached to the centre mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall. The integrated gas cylinders enable the bench to be folded out of the way when not in use to maximise floor space in front of the unit.</p> <p>Electric with wired hand control – the red buttons clearly indicate “Up” and “Down”. The wired hand control comes with a holder for fitting to the wall.</p> <p>The stretcher of the bench is very easy to clean, and the stretcher cover can be unzipped for washing or easy replacement.</p>
Mains connection (V)	120
Actuator (motor)	LA31 Careline, (Linak), 6000 N, IPX6
Control box	CB9 Careline, (Linak), IPX6
Wired remote control	HB3 (Linak), IP66
Power consumption (maximum load)	<p>Power consumption by maximum load, watt: 70 W</p> <p>Power consumption by standby: 2 W</p>
Warranty	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.</p>



Product	Nursing bench 2000, coated mesh cover, electrically height adjustable
Item number	R8468000-93
Cleaning	<p>SCT 1000 & 2000</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes OxiGenesis from ENVIROX</p> <p>NON-DISINFECTANTS:</p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>



PRESSALIT®

Product	Nursing bench 2000, coated mesh cover, electrically height adjustable
Item number	R8468000-93

Installation in different wall types

	<p>Concrete >B25</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8658 - Mounting kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.</p>	



PRESSALIT®

Product Item number	Nursing bench 2000, coated mesh cover, electrically height adjustable R8468000-93
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.</p>
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm)</p> <p>Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.</p>
	<p>Drywall (2 x 12,5 mm) water-resistant</p> <p>Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.</p>

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®