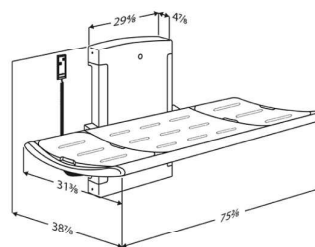




PRESSALIT®

Product Item number	Adult changing table 3000, electrically height adjustable R8520021
Color	EAN: 5708590265629 FDA: INQ/D124319 Lead grey
Product description	Discontinued per 31 Dec 2018 and replaced by R8539 Adult changing table 3000, powered, height adjustable 27.6"; from 11.8" - 39.4", with wired hand control Foldable. Incl. water collection tray with integrated outlet. Surface length 70.9", outer measurements 38.9" x 75.2". Excl. mounting kit. See Technical specifications (pdf) for information on mounting kits suitable for various wall types. Discontinued per 31 Dec 2018 and replaced by R8539 Adult changing table 3000, powered, height adjustable 27.6"; from 11.8" - 39.4", with wired hand control Foldable. Incl. water collection tray with integrated outlet. Surface length 70.9", outer measurements 38.9" x 75.2". Excl. mounting kit. See Technical specifications (pdf) for information on mounting kits suitable for various wall types.
Adjustability	Height adjustable
Adjustment height	27,58 inch
Product capacity	Maximum load: 440 lbs. When the table is loaded with 440 lbs, each mounting bolt supports 57 lbs. <i>This product has been tested statically with a load of 1100 lbs as the IEC 60601-standard requires testing with a safety factor.</i>
Speed up	Speed up: 0,59 inch/s Speed down: 0,79 inch/s
Noise level	48 dB
Materials	Aluminium, polyurethane foam, plastics, powder-coated steel





Product Item number	Adult changing table 3000, electrically height adjustable R8520021
Ergonomy Features	Pressalit offers a range of high quality, fixed or power adjusted changing tables designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centers, institutions, hospitals and private homes. The changing table 3000 series provides maximum safety and comfort. Besides an award winning design, the series includes comfortable and durable surfaces, including polyurethane foam and aluminum. The flexibility of the changing table provides optimal working conditions for the caregiver as well as a high level of comfort and safety for the user. A water collection tray with integrated water outlet prevents water from splashing on caregiver or the floor. The flexible hose attached to the center mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall. Adjustable neck support which can be moved to either end of the 70.9" changing tables. The slightly curved form makes the table extremely comfortable and safe. The integrated pneumatic cylinders enable the table to be folded out of the way when not in use to maximize floor space.
Mains connection (V)	120
Actuator (motor)	LA31 Careline, (Linak), 6000 N, IPX6
Control box	CB9 Careline, (Linak), IPX6
Wired remote control	HB3 (Linak), IP66
Power consumption (maximum load)	Power consumption by maximum load: 70 W Power consumption by standby: 2 W
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.



PRESSALIT®

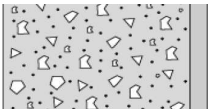
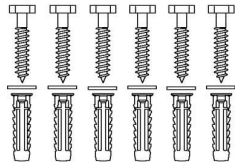
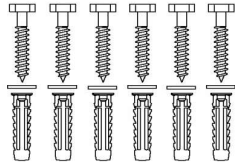
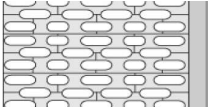
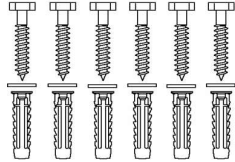
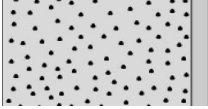
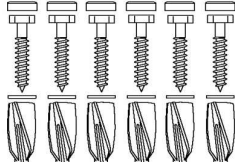
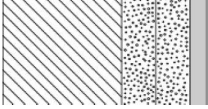
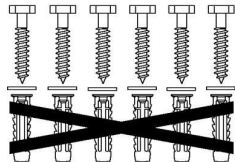
Product Item number	Adult changing table 3000, electrically height adjustable R8520021
Cleaning	This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes OxiGenesis from ENVIROX NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate. DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.



PRESSALIT®

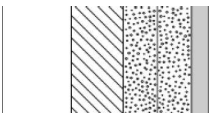
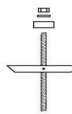
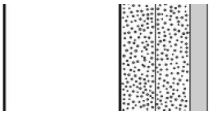
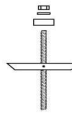
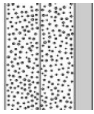
Product	Adult changing table 3000, electrically height adjustable
Item number	R8520021

Mounting in different wall types

	Concrete >B25 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Brick >= Mz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.	
	Aerated Concrete P4 V8658 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs GB14. To be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.	



PRESSALIT®

Product Item number	Adult changing table 3000, electrically height adjustable R8520021		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com